
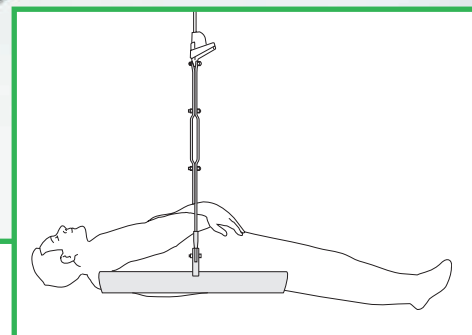


# Rigor

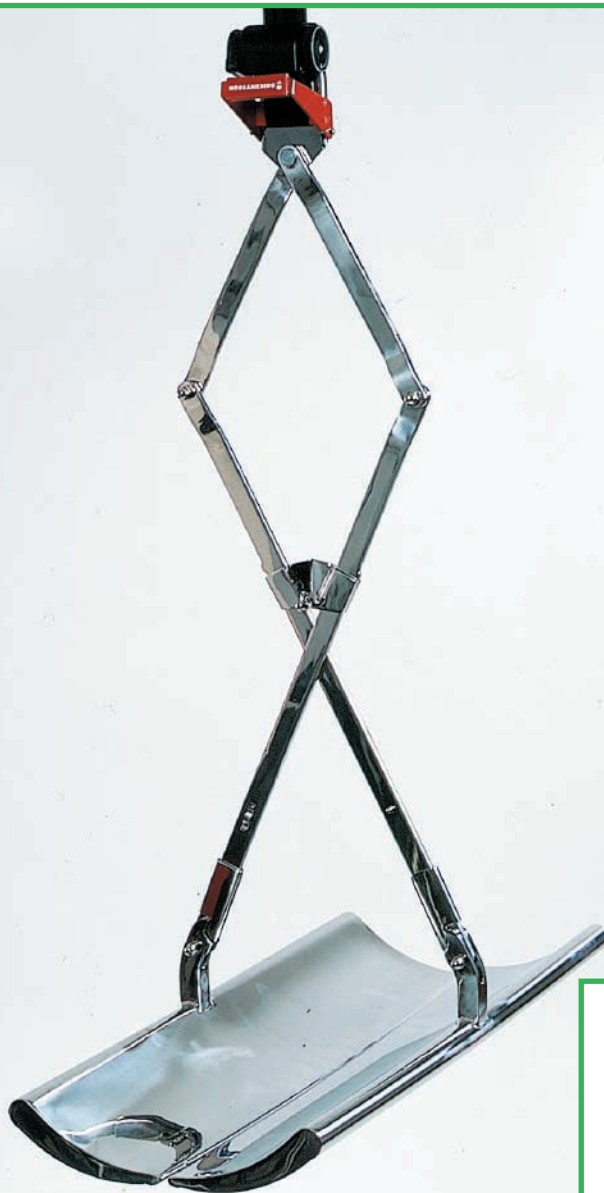

Customized lifting assistance for a special environment. With Rigor, one person can easily, conveniently and safely transfer a body, even where space is limited.

The scissor-like design, entirely in stainless steel, ensures maximum safety. The lift-blades are specially designed to fit securely around the body. In an elevated position, they are locked individually, but the locking device is easily opened so that the lift-blades can be released when the body is placed in a coffin. Rigor is hygienic, easy to clean, and fits all of Liko's stationary lifts.

Weight: 10.3 kg (22.66 lbs)    Width: 35 cm (13¾")    Length: 75 cm (29½")    Total height: 96 cm (37¾")

Prod. No. 3156026



Maximum load: 150 kg (330 lbs.)

Made in Sweden.

Liko is quality certified according to ISO 9002 and environment certified according to ISO 14001.



 **Liko**

[www.liko.com](http://www.liko.com)

Manufactured by  
**Liko AB**  
Nedre vägen 100  
SE-975 92 Luleå, Sweden  
info@liko.se

Liko Distributor in UK  
**Liko (UK) Ltd.**  
Liko House, L&M Business Park  
Norman Road, Altrincham Cheshire  
WA 14 4ES, UK  
Tel (0161)926 8067 Fax (0161)929 7050  
info@liko.co.uk

National Distributor in USA and Canada  
**Liko Inc**  
842 Upper Union Street, Suite 4  
Franklin, Ma 02038, USA  
Tel 508-553-3993, 1-888-545-6671  
Fax 508-528-6642  
info@likoinc.com