

Ultra Twin FreeSpan



Ultra Twin on a FreeSpan railsystem – a freestanding overhead lift mounted on a multiple caster chassis. With a lifting capacity of 400 kg (880 lbs.), Ultra Twin FreeSpan is designed for the heaviest lifting situations. The two interacting Likoall 242 ES R2R lift motors make it possible for the assistant to both lift and position the patient while in the sling during lift and transfer. The two lift motors may be simultaneously operated with one hand control, or used with separate hand controls with the stretcher unit.

Ultra Twin FreeSpan



Ultra Twin FreeSpan comes in two different models: For lifts with a sling or for horizontal lifts.

Ultra Twin FreeSpan is available with different rail lengths, ranging from 2.0-5.5 meters (6' 7½"–18'). Choose the length depending on needs and the space available. UltraTwin FreeSpan is easy to assemble and move and the rails can be raised or lowered when needed.


Prod.No.

Products needed for a complete Ultra Twin sling unit:

3103520	Ultra Twin side supports	1 pair
3103115-50	FreeSpan Liftrail LR	2 pcs
3122007	Likorall 242 ES R2R	2 pcs
3126040	Ultra control unit	1 set
3126101-4	Battery charger (EU, UK US, AU)	1 pce
3103603	Ultra Slingbar	1 pce
	Appropriate sling/accessories	

Products needed for a complete Ultra Twin stretcher unit:

3103521	Ultra Twin sidesupport, wide	1 pair
3103115-50	FreeSpan Liftrail LR	2 pcs
3122007	Likorall 242 ES R2R	2 pcs
3126034	Handcontrol LR 242 ES, 2 button	2 pcs
3126101-4	Battery charger (EU, UK US, AU)	2 pcs
3156058	Ultra Stretcher	1 pce
	Appropriate sling/accessories	

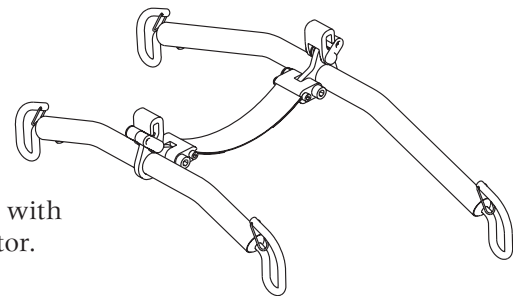

Ultra Slingbar with safety connector.

Prod. No.

3103603	Ultra Slingbar
---------	----------------



Maximum load: 400 kg, 880 lbs.

Made in Sweden

Liko is quality certified according to ISO 9002 and environment certified according to ISO 14001.



www.liko.com

Manufactured by
Liko AB
Nedre vägen 100
SE-975 92 Luleå, Sweden
info@liko.se

Liko Distributor in UK
Liko (UK) Ltd.
Liko House, L&M Business Park
Norman Road, Altrincham Cheshire
WA 14 4ES, UK
Tel (0161)926 8067 Fax (0161)929 7050
info@liko.co.uk

National Distributor in USA and Canada
Liko Inc
842 Upper Union Street, Suite 4
Franklin, Ma 02038, USA
Tel 508-553-3993, 1-888-545-6671
Fax 508-528-6642
info@likoinc.com